
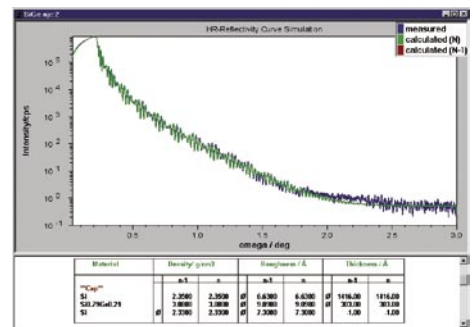


Changer for Attenuation Filters

High Speed Automatic Changer of Attenuation Filters for high intensity dynamics

Features

- Count rates from below 0.1 cps to above 10^9 cps
- For use in XRD 3003 Diffraction systems (mounted on adjustment base)
- Stepper motor driven, controlled by 1 axis of C 3006 controller
- Fully automatic control with Rayflex Software

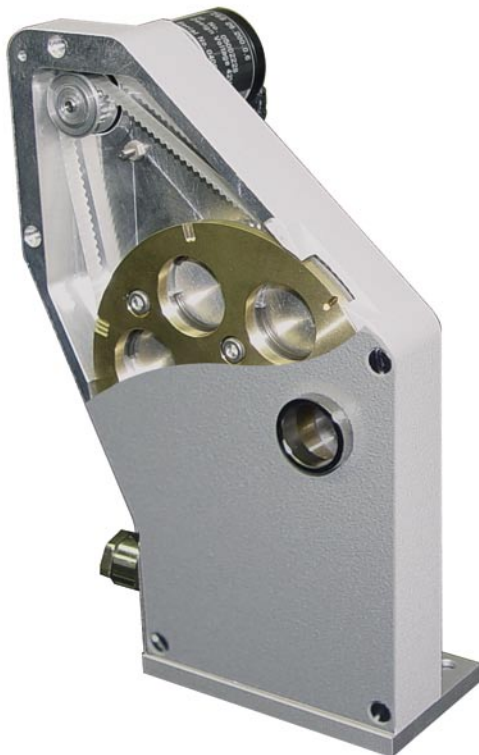

Example from a Reflectometry Measurement

6 Filter Wheel Positions:

- Pos. 1 – open
- Pos. 2 - 0.15 mm Al
- Pos. 3 - 0.30 mm Al
- Pos. 4 - 0.45 mm Al
- Pos. 5 - 0.60 mm Al
- Pos. 6 - 0.75 mm Al

